



Disposable Veress Needles


Fairmont Medical's Veress Needles are designed for creating pneumoperitoneum and feature a sharp tip for easy insertion.

Fairmont Medical's, Disposable Veress Needles are designed to create a pneumoperitoneum in laparoscopic procedures. Our unique red indicator window design incorporates a retracting spring loaded blunt stylet and sharp hollow bore needle.

Used to establish pneumoperitoneum prior to trocar and cannula insertion, our Veress Needle is available in 120mm, 150mm and 200mm lengths.

- Safe blunt stylet shields needle tip once correct positioning has been achieved
- Easy to view windowed indicator and spring-loaded stylet provides sensory feedback for correct needle placement
- Sharp beveled needle tip ensures simple, safe abdominal wall penetration
- High-flow stopcock with luer lock connector for convenient connection of insufflation tubing
- Easy to identify colour coding
- Single use, sterile, ready to use
- Quality Australian designed and manufactured

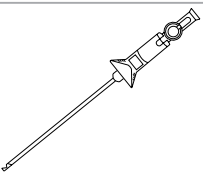

Sterile, Disposable				
Product	Product Code	Colour Grip	Length	Box Qty
	DVN7012	Blue	120mm	30
	DVN7015	Yellow	150mm	30
	DVN7200	Orange	200mm	30



Corporate Office
 Fairmont Medical
 11 Scoresby Road Bayswater 3153 VIC
 Ph 1300 972 088
 Ph International +613 9720 8840
 Fx International +613 9237 9399

UK Head Office
 Fairmont Medical Products Ltd
 Elizabeth House
 28 Baddow Road Chelmsford
 Essex CM2 0DG
 Ph: +44 (0)1245 206800

