



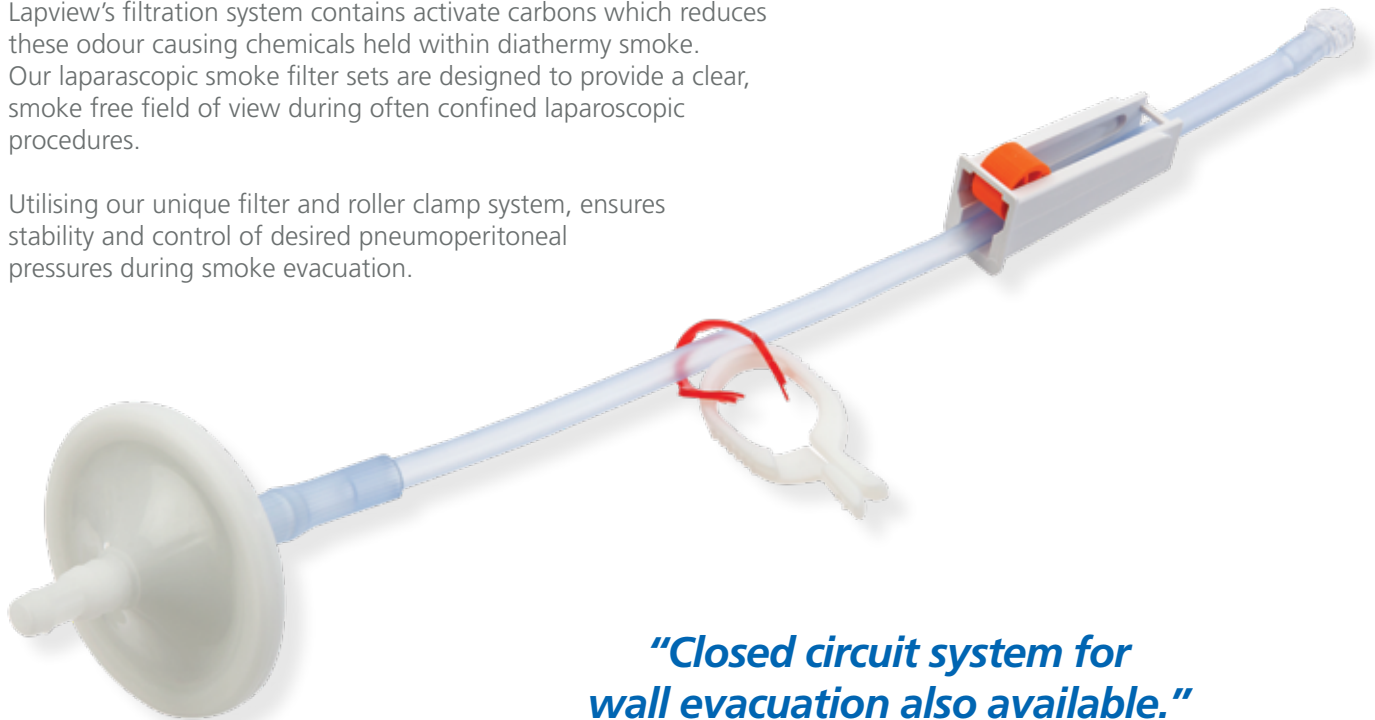
# Laparoscopic Smoke Filter

Fairmont Medical's Lapview is a passive and on demand smoke evacuation filter system which provides a safe, cost effective smoke removal solution.

For decades Fairmont Medical has been at the forefront of providing a safe and efficient method of removing diathermy smoke during minimal invasive surgery. Smoke created during the burning of tissue with diathermy instruments, contains chemicals which are known irritants and in the case of PAH's, benzene and formaldehyde, known carcinogens.

Lapview's filtration system contains activate carbons which reduces these odour causing chemicals held within diathermy smoke. Our laparoscopic smoke filter sets are designed to provide a clear, smoke free field of view during often confined laparoscopic procedures.

Utilising our unique filter and roller clamp system, ensures stability and control of desired pneumoperitoneal pressures during smoke evacuation.



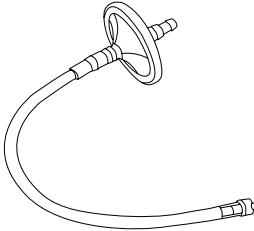
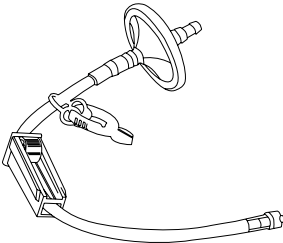
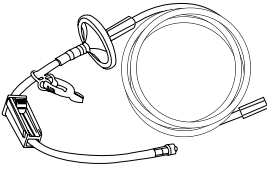
***"Closed circuit system for wall evacuation also available."***

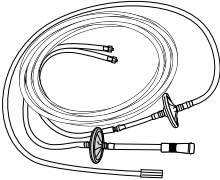
- Safe smoke removal for both patient and theatre personnel
- Reduces patient risk of mutagenic chemical absorption
- Optional controlled passive or active diathermy smoke removal maximum 4L/min flow rate
- Protects theatre personnel against potentially infectious blood-borne pathogens utilising our bacterial viral PTFE filter media
- Easy to connect ISO standard luer lock system to trocar ports
- Saving valuable operating time
- Closed circuit system also available
- Single use, sterile, ready to use
- Phthalate (DEHP) and latex free
- Quality Australian designed and manufactured

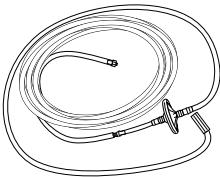




# Laparoscopic Smoke Filter

Lapview - Laparoscopic Smoke Filter Set - Sterile, Disposable						
Product	Product Code	Filters	Tubing ID	Tubing OD	Length	Box Qty
	LESS7018	2x Activated Carbon Filter 59gsm PTFE Filter (Maximum 4 litre extraction rate)	5.4mm	7.2mm	300mm	40
	LESS7008	2x Activated Carbon Filter 59gsm PTFE Filter with roller clamp	5.4mm	7.2mm	300mm	40
	LESS7006	2x Activated Carbon Filter 59gsm PTFE Filter with roller clamp	5.4mm	7.2mm	600mm	40
	LESS7007	2x Activated Carbon Filter 59gsm PTFE Filter with roller clamp	Filter Tubing 5.4mm Suction Tubing 6mm	Filter Tubing 7.2mm Suction Tubing 9mm	Filter Tubing 600mm Suction Tubing 3m	40

Dual CO <sub>2</sub> Filter And Tubing / In-Line Smoke Filter Set - Sterile, Disposable						
Product	Product Code	Filters	Tubing ID	Tubing OD	Length	Box Qty
	LESS7025	CO <sub>2</sub> Filter & Smoke Filter	6mm	9mm	3m	25

Laparoscopic Evacuation Smoke Set - Lapview Smoke/Plume Filter With 6 m Tubing Luer/Hex Connector - Sterile, Disposable						
Product	Product Code	Filters	Tubing ID	Tubing OD	Length	Box Qty
	LESS7030	Smoke Filter	6mm	9mm	6m	35



Corporate Office  
Fairmont Medical  
11 Scoresby Road Bayswater 3153 VIC  
Ph 1300 972 088  
Ph International +613 9720 8840  
Fx International +613 9237 9399

UK Head Office  
Fairmont Medical Products Ltd  
Elizabeth House  
28 Baddow Road Chelmsford  
Essex CM2 0DG  
Ph: +44 (0)1245 206800

